



Intended use:

AflaStar™ R immunoaffinity columns contain monoclonal antibodies against aflatoxin B1, B2, G1 and G2, which are covalently bound to gel-particles and are intended for the analysis of various commodities (e.g. cottonseeds, cereals, peanuts, walnuts, pistachios, spices).

Recommended solvents and buffers:

- Extraction solution: 84/16 (v/v) acetonitrile/water or 60/40 (v/v) acetonitrile/water or 60/40 (v/v) methanol/water (all HPLC grade)
- Dilution: PBS-buffer (10 mM)

PBS-buffer:

- 8 g NaCl, 1.2 g Na_2HPO_4 , 0.2 g KH_2PO_4 , 0.2 g KCl (p.A.)
- dissolve in 990 mL distilled or deionized water
- adjust pH to 7.4 using NaOH (1 M) or HCl (1 M)
- fill up to 1000 mL with distilled or deionized water
- Rinse solution: distilled water or deionized water or PBS-buffer (10 mM)
- Eluent: methanol (HPLC grade)

All solvents and buffers should be at room temperature (15 – 25 °C).

Romer Labs recommends the use of Biopure™ isotope labeled internal standards.

Storage:

Always store at 2 – 8 °C (35 – 46 °F) when not in use.

Do not freeze. Do not use the IAC beyond the expiration date.


IMPORTANT: Download certificate of analysis by scanning the QR code on the external label or by visiting <https://www.romerlabs.com/en/customer-resources/>

Special Notes for column use:

- StarLine™ IAC contain sodium azide.
- StarLine™ IAC are designed for single use only.
- Application notes for the use of StarLine™ IAC have been developed by Romer Labs and are available on request.

Disclaimer/Warranty:

The user or buyer assumes all risk in using Romer Labs products and services. Romer Labs makes no warranty of any kind, either expressed or implied, except that the materials from which its products are made are of standard quality. If any materials are defective, Romer Labs will, at its own discretion and option, repair or replace any product or components, or repeat services, which are proven to be defective in workmanship or material within its warranty periods. There is no warranty of merchantability of the products, or of the fitness of the product for any purpose. Romer Labs shall not be liable for any damages, including special or consequential damage, lost profits or expenses arising directly or indirectly from the use of its products or services. Excepting the aforementioned, Romer Labs hereby disclaims all other remedies, warranties, guarantees or liabilities, expressed or implied, arising by law or otherwise.

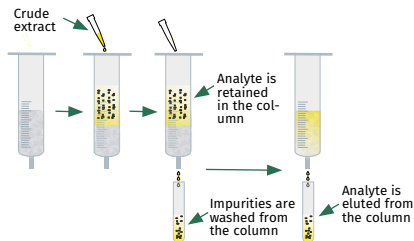

Procedure

1 Extraction

- Weigh 25 g of a representative ground sample (0.5 mm) into a suitable container.
- Add 100 mL of extraction solution (e.g. 86/16 (v/v) acetonitrile/water).
- Blend at high speed for 3 minutes or shake for 1 hour on an orbital or gyratory shaker.
- Centrifuge or filter the supernatant through qualitative filter paper.
- Liquid samples can be diluted directly with PBS or mixed with solvent to obtain the appropriate ratios of the reagent

2 Dilution

- used (e.g. if using 84/16 AcN/H₂O (v/v): 0.640 mL of sample and 3.36 mL of AcN).
- Any solvent content in the sample must be considered in the calculation of the dilution.
- Dilute the extract with PBS until the content of acetonitrile is not higher than 10% (v/v). For MeOH extracts the content of MeOH can be up to 20% (v/v).
- Check pH and adjust with NaOH to 6 – 8.

3 Sample Application

- The IAC must be at room temperature (15 – 25 °C) for usage!
- It is not necessary to rinse the IAC before applying the extract!
- Apply the diluted extract to the IAC.
- Let all extract pass through the IAC with a flow rate of approx. 1 – 3 mL/min.

4 Rinse

- Rinse the IAC with 2x10 mL distilled water, deionized water or PBS at a flow rate of 1 – 3 mL/min.

5 Elution

- Place a suitable vial under the AflaStar™ R IAC.
- Use 1.5 – 3.0 mL of methanol for the elution of bound aflatoxins; the eluent should be applied to the column in several small portions (e.g. 3x0.5 mL)!
- Leave the eluent on the column for a few seconds before starting elution to allow intensive contact with the gel.
- Remove any remaining liquid from the IAC by applying pressure at the top or vacuum at the bottom.
- In case of low level contamination the eluent can be dried down and re-dissolved in a small portion of mobile phase (use silanized glasware).
- Inject

You can find worldwide contact information and learn more about our complete line of products for mycotoxin testing on our website.