

Sterile devices for sampling



The sampling materials

- 1 A sampling material: Swab cloth or SteriSox®
- 2 Premoistened or dry
- 3 Packed in a sampling bag which can be closed
- 4 Available with or without accessories




The sampling bags


1 Sampling materials for surfaces

Sodibox® Swab cloth

- Swab cloth size > 900 cm² (cf. EU regulation)
- Does not drip, fray or tear
- Unwoven, consistent smooth or honeycombed structure
- No bactericidal treatment

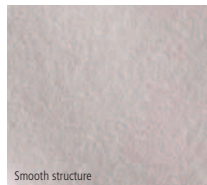

Smooth structure

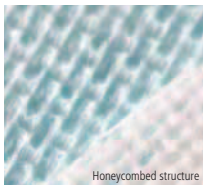

Honeycombed structure

SteriSox®

- Knitted jersey, stretchable and strong
- Easy to put on with wide opening
- Elasticated top



Knitted jersey

You can put up to 5 SteriSox® over a size 10(UK), 11(US), 45(FR) boot with over-boots!



Effectiveness:

- Sample is representative of the tested surface
- Collects quality sample with strong fixation of germs
- Fast release of sample in pre-enrichment broth



2 Diluents

Distilled water or Tryptone Salt

For testing during production.

Premoistened devices are obligatory for the regulated sectors of Gallus gallus and turkey.

For the other sectors such as duck, Guinea fowl, pig, we produce devices without diluent (dry).

Buffered peptone water at 10% neutralizing agent

For testing after disinfection.

Reliability:

- Reproducible and calibrated dose
- Volume proportional to the size and absorption level of the sampling material
- Diluents made in our laboratory (XP CEN ISO/TS 11133)

3 Sampling bag

DESIGNED FOR EASY USE

- Easy opening with pre-cut top
- Integral closing system
- Large enough to collect the necessary amount of pre-enrichment solution
- Wide with resistant seams for use in a stomacher®



4 Packaging

Single unit: with accessories

Device with accessories in a zip bag.

Advantages: Complete kit. Samples can be packed individually, limiting the risk of cross contamination.




SteriSox® + 1 pair of gloves and 1 pair of over-boots



Swab cloth + 1 pair of blue vinyl gloves



SteriSox® economical

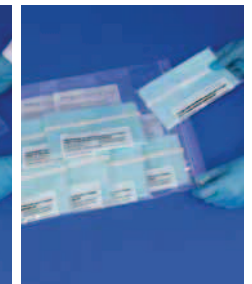

Swab cloth economical

Economical: without accessories

Total of 15 (ECO 15) or 25 (ECO 25) devices in a outer zip bag. See page 7 "Storage advice".

Accessories for Economical packaging



Ref. 4050

Ref. 4050 • 2 gloves and 2 sterile over-boots for SteriSox®
Packaged in a zip bag 18 x 25 cm, over-boots with elasticated top, and blue vinyl gloves, size 8-9(Medium), box of 100 units.

Ref. 4004 • 2 sterile gloves for Swab cloth
Packaged in a zip bag 18 x 25 cm, Blue vinyl gloves, size 8-9(Medium), box of 350 units.



Ref. 4004

Device for hatcher basket liners (cf. regulation EU n° 213/2009 and n° 200/2010)

Single unit packaging:

- Zip bag dim. 300 x 400 mm with batch number, containing:
 - One stomacher® type bag dimensions 380 x 510 mm, with identification label, that will receive at least 5 soiled hatcher basket liners.
 - One pair of blue vinyl gloves size 8-9 (Medium) Box of 200 units.



Ref. 4003 - Device for hatcher basket liners

